# Document made available under the Patent Cooperation Treaty (PCT)

International application number: PCT/EP2006/002172

International filing date: 09 March 2006 (09.03.2006)

Document type: Certified copy of priority document

Document details: Country/Office: EP

Number: 05005390.9

Filing date: 11 March 2005 (11.03.2005)

Date of receipt at the International Bureau: 20 June 2006 (20.06.2006)

Remark: Priority document submitted or transmitted to the International Bureau in

compliance with Rule 17.1(a) or (b)





Europäisches **Patentamt** 

European **Patent Office**  Office européen des brevets

Bescheinigung

Certificate

Attestation

Die angehefteten Unterlagen stimmen mit der ursprünglich eingereichten Fassung der auf dem nächsten Blatt bezeichneten europäischen Patentanmeldung überein.

The attached documents are exact copies of the European patent application of photomes a la version described on the following sinitialement déposée de page, as originally filed.

la demande de brevet jeuropéen spécifiée à la page suivante.

Patentanmeldung Nr.

Patent application No. Demande de brevet nº

05005390.9

Der Präsident des Europäischen Patentamts; Im Auftrag

For the President of the European Patent Office

Le Président de l'Office européen des brevets p.o.

R C van Dijk



Europäisches Patentamt European Patent Office Office européen des brevets

PCT/EP200 6/ 0 0 2 1 7 2

Anmeldung Nr:

Application no.: 05005390.9

Demande no:

Anmeldetag:

Date of filing:

Date de dépôt:

11.03.05

Anmelder/Applicant(s)/Demandeur(s):

AXXAM S.r.l. Via Massena, 12/7 20145 Milano ITALIE

Bezeichnung der Erfindung/Title of the invention/Titre de l'invention: (Falls die Bezeichnung der Erfindung nicht angegeben ist, siehe Beschreibung. If no title is shown please refer to the description. Si aucun titre n'est indiqué se referer à la description.)

Photoproteins with enhanced bioluminescence and assays using the same

In Anspruch genommene Prioriät(en) / Priority(ies) claimed /Priorité(s) revendiquée(s)
Staat/Tag/Aktenzeichen/State/Date/File no./Pays/Date/Numéro de dépôt:

Internationale Patentklassifikation/International Patent Classification/Classification internationale des brevets:

C07K14/00

Am Anmeldetag benannte Vertragstaaten/Contracting states designated at date of filing/Etats contractants désignées lors du dépôt:

AT BE BG CH CY CZ DE DK EE ES FI FR GB GR HU IE IS IT LT LU MC NL PL PT RO SE SI SK TR LI

10

15

20

25

## PHOTOPROTEINS WITH ENHANCED BIOLUMINESCENCE AND ASSAYS USING THE SAME

The present invention relates to photoproteins with enhanced bioluminescence and to their use as intracellular calcium indicators. The photoproteins are obtained by mutagenesis of the clytin coding sequence and show enhanced bioluminescence, high affinity for calcium and long-lasting light emission. They are conveniently used in cell-based assays to determine variations of intracellular calcium concentration, particularly in assays for the screening of molecules with high and ultra-high-throughput techniques.

## BACKGROUND OF THE INVENTION

Bioluminescence is the phenomenon by which visible light is emitted by living organisms or by a substance derived from them through a variety of chemiluminescent reaction systems. Bioluminescence reactions require three major components: a luciferin, a luciferase and molecular oxygen. However, other components may also be required in some reactions, including cations (Ca++ and Mg++) and cofactors (ATP, NAD(P)H). Luciferases are enzymes that catalyse the oxidation of a substrate, luciferin, and produce an unstable intermediate. Light is emitted when the unstable intermediate decays to its ground state, generating oxyluciferin. There are many different unrelated types of luciferin, although many species from at least seven phyla use the same luciferin, known as coelenterazine. In some animals (e.g. jellyfish) the luciferin/luciferase system can be extracted in the form of a stable "photoprotein" which emits light upon calcium binding. Photoproteins differ from luciferases in that they are stabilized oxygenated intermediate complexes of luciferase and luciferin. Photoproteins are present in many marine coelenterates and allow these organisms to emit light for a variety of purposes including breeding, feeding and defense (1). There are many luminescent

10

15

20

25

organisms, but only seven photoproteins, namely Thalassicolin (2,3), Aequorin (4-6), Mitrocomin (syn. with Halistaurin) (7,8), Clytin (syn. with Phialidin) (8,9), Obelin (2,6,10,11), Mnemiopsin (12,13) and Berovin (12,13) have been isolated so far. All these proteins are complexes formed by an apoprotein, an imidazopyrazine chromophore (coelenterazine) and oxygen. Their structures are highly conserved, especially in the region containing the three calcium binding sites (EF-hand structures). The term "photoprotein" identifies the luciferin-bound polypeptide, which is capable of luminescence, while "apophotoprotein" is used to indicate the protein without luciferin.

The most studied photoproteins are Aequorin, isolated from Aequorea victoria (14) and Obelin, isolated from Obelia longissima (15). The photoprotein may be regenerated from the apophotoprotein by incubation with coelenterazine, molecular oxygen, EDTA and 2-mercaptoethanol or dithiothreitol. Since coelenterazine is the common luminescent substrate used by the photoproteins Aequorin, Mitrocomin, Clytin and Obelin, the light-emitting reaction is likely the same in these four photoproteins (16,17).

The Clytin photoprotein was cloned in 1993 by Inouye et al. (18). To date not much work has been done on this photoprotein. The primary structures of aequorin, mitrocomin, clytin and obelin were aligned and showed very strong amino acid sequence identities. The Ca<sup>2+</sup>-binding sites of Clytin were also found to be highly conserved (19). It was found that hydrozoan Ca<sup>2+</sup>-binding photoprotein differs from other Ca<sup>2+</sup>-binding proteins such as calmodulin and troponin C by a relatively high content of cysteine, histidine, tryptophan, proline and tyrosine residues.

The analysis of the primary structure of clytin shows that it contains 198 aminoacidic residues (aa) and belongs to the family of photoproteins.

Photoproteins are widely used in reporter gene technology to monitor the cellular events associated with signal transduction and gene expression.

10

15

20

25

The study of cellular events and their regulation requires sensitive, non invasive analytic methods. Photoproteins and in general the use of bioluminescence are excellent reporter systems as they have virtually no background in contrast to fluorescence systems.

Photoproteins are expressed in mammalian cells to monitor calcium changes in response to different stimuli. Intracellular calcium concentrations can be measured by adding the cofactor coelenterazine to mammalian cells expressing the photoprotein and detecting photon emission, which is indicative of intracellular calcium concentration. The use of cells which express both a photoprotein and a receptor involved in the modulation of intracellular calcium concentration provides a valid system for the screening of compounds for their effects on the release of intracellular calcium. High throughput screening assays can also be designed using a photoprotein as reporter system. The sensitivity of the system as well as its high signal to noise ratio allow the use of small assay-volumes. Aequorin is up to now the most used photoprotein for these screening assays.

Calcium flux assays are commonly carried out in HTS format utilizing optical screening apparatuses suited for the simultaneous analysis of a high number of samples and equipped with a luminescence imaging system with a CCD Camera detector. However, one of the most used instruments in HTS is FLIPR® (Fluorometric Imaging Plate Reader, Molecular Devices Corporation, Sunnyvale, CA, USA) which was developed as a high throughput optical screening tool for cell-based fluorescent assays. The apparatus is equipped with an optical detection device that allows for signal isolation on a cell-monolayer, thereby enhancing sensitivity for cell-based assays. The excitation source can be either an Argon laser or a broadband source as a Xenon lamp.

With a light-tight enclosure, extremely sensitive and fast camera, and

true simultaneous on-line liquid dispensing, the FLIPR® system most recent versions (FLIPR³ and FLIPRTETRA) have been made suitable also for luminescence assays, even if with lower sensitivity compared to CCD Camera-based equipments.

For the use of the above described instruments FLIPR<sup>®</sup>, FLIPR<sup>3</sup> and FLIPR<sup>TETRA</sup> and in general for all instruments with a low sensitivity for luminescence assays, a photoprotein with enhanced light emission is highly desirable.

#### DISCLOSURE OF THE INVENTION

According to a first aspect, the invention provides an isolated photoprotein containing an amino acid sequence which:

- a) is able to bind coelenterazine and calcium, producing bioluminescence;
- b) has an identity of at least 90%, preferably of at least 95%, more preferably of at least 98% to SEQ ID NO: 1 (Clytin);
- c) in sequence alignment with SEQ ID 1 (Clytin), presents one of the following single or multiple substitutions (the residue positions are referred to SEQ ID NO: 1):
- i)  $C_{54} \rightarrow S$ ;
- 20 ii)  $S_{132} \rightarrow C$ ;

5

15

- iii)  $K_{48}\rightarrow R$ ,  $N_{195}\rightarrow D$ ;
- iv)  $Q_{68}\rightarrow R$ ,  $A_{90}\rightarrow V$ ,  $T_{184}\rightarrow I$ ;
- v)  $Y_{82} \rightarrow F$ ,  $K_{110} \rightarrow N$ ,  $F_{125} \rightarrow L$ ,  $S_{149} \rightarrow R$ ;
- vi)  $G_{142} \rightarrow C$ ;
- 25 vii)  $I_{53} \rightarrow V$ ,  $S_{149} \rightarrow R$ ;
  - viii)  $N_{18}\rightarrow D$ ,  $I_{40}\rightarrow V$ ,  $K_{56}\rightarrow R$

In a preferred embodiment, the photoprotein contains an amino acid sequence which is selected from the group of SEQ ID NO: 2, 3, 4, 5, 6, 7, 8

10

15

20

and 9. Compared to known or commercially available photoproteins, the photoproteins of the invention show improved bioluminescence activity, higher affinity to calcium and longer-lasting light emission.

Apart from the indicated residue-substitutions, which confer the desired photoprotein bioluminescence activity, the Clytin sequence can be further modified without negatively affecting the photoprotein's bioluminescence activity, especially by conservative substitutions of amino acidic residues, within the indicated sequence-identity limits. In addition, the Clytin sequence can be deleted of small portions, without altering its photoprotein activity.

In a further aspect, the invention is directed to a polynucleotide encoding a photoprotein as defined above. In a preferred embodiment, the polynucleotide sequences are optimized for mammalian codon usage according to SEQ ID NO: 10, 11, 12, 13, 14, 15, 16, 17. In a further preferred embodiment, the nucleic acid molecules are fused to mitochondrial target sequences (20, 21, 22).

According to a further aspect, the invention provides expression vectors and host cells containing the indicated polynucleotides. Host cells expressing a photoprotein according to the invention produce an intense bioluminescence in response to calcium stimulation, which is much higher than that observed with natural photoproteins, in particular with the most used one, Aequorin.

In a further aspect, the invention provides a cell-based assay for determining intracellular calcium concentration by means of a photoprotein according to the invention.

In a preferred embodiment, the changes in intracellular calcium concentration are determined by:

- a) providing a cell expressing a photoprotein of SEQ ID NO: 2-9, variants or fragments thereof;
- b) loading of the cells with coelenterazine;

10

15

20

25

- c) contacting the cells with an agent stimulating calcium influx or calcium release from intracellular stores;
- d) detecting the photoprotein's bioluminescence.

The assay is preferably carried out in a high-throughput format utilizing an optical screening tool or apparatus suited for multi-sample analysis, such as a luminescence imaging system with a CCD Camera detector for high and ultra high throughput applications, or with the Fluorometric Imaging Plate Reader (FLIPR\*). With both these systems, the photoproteins of the invention produce the highest signal compared to known photoproteins commonly used in automatized cell functional assays.

In a preferred embodiment, cells expressing a photoprotein and a receptor involved in intracellular calcium mobilization are used to test candidate molecules for their effects on receptor modulation. Typically, cells are transfected with an expression vector containing a photoprotein encoding sequence and when not endogenously present, a receptor or channel of interest. The positive clones are selected and plated in a suitable medium, cultured cells are loaded with the coelenterazine substrate and the assay is started by adding the test molecule or stimulus. The produced luminescence is read by a suitable detection system (CCD camera or luminometer). The assay can also be run in an automatic apparatus equipped with multi-well plate reading, in particular the FLIPR® system. In this case, photoprotein-expressing cells are plated in microplate wells, which, after addition of the test molecule/stimulus, are read simultaneously with signal recording.

High throughput screening assays set up with a photoprotein-based reporter system show improved sensitivity and signal-to-noise ratio. Cells expressing a photoprotein of the invention produce an intense bioluminescence in response to calcium stimulation, which is generally higher than that observed with natural photoproteins.

In a further aspect, the invention provides an assay kit containing a preparation of cells expressing an invention photoprotein under the control of a stable or inducible promoter, and reagents suitable for running the assay.

In addition, the photoproteins of the invention may be used as intracellular calcium indicators in diagnostic methods based on the measurement of cellular calcium ion concentration and/or cellular calcium ion influx/outflow.

The invention will be described in more detail in the following experimental section.

#### MATERIALS AND METHODS

#### Reagents

5

10

15

20

25

Restriction enzymes were purchased from New England Biolabs and used according to supplier's instructions. The Ligation Independent Cloning (LIC) kit was from Novagen (Nottingham, UK). For in vitro Transcription and Translation we used the TNT Quick coupled kit from Promega (Madison, WI). Reagents for PCR, and competent cells of E. coli strains XL-1Blue and BL21-Gold(DE3), were from Stratagene (La Jolla, CA). Oligonucleotides were purchased from Primm (Milan). Coelenterazine was from Pharma Tech. International Inc. (Fairfield, NJ). All other chemicals were from standard sources and were of reagent grade or better.

- 1. Generation of a Randomly Mutagenized Library and Screening
- 1.1 Photoprotein optimization for expression in mammalian cells (GENEART GmbH, Regensburg, Germany)

The codon usage of the wild-type clytin gene was adapted to the codon bias of highly expressed mammalian genes. In addition regions of very high (> 80%) or very low (< 30%) GC content have been avoided where possible.

For efficient translation initiation the Kozak-consensus sequence was introduced upstream of the start codon. Two STOP codons were added to

ensure efficient termination.

5

10

25

#### 1.2 Random Mutagenesis

The GeneMorph II Random Mutagenesis kit (Stratagene) was used following supplier's instructions. Two different initial amounts of target DNA were used to achieve a high mutation rate, 0.lng and 0.0lng.

PCR primers were appropriately designed to contain 5' LIC extensions (in italics) corresponding to sequences described in the Ek/LIC Cloning Kit (Novagen)

Upper: GATGACGACGACAAG-ATGGCCGACACCGCCAG

Lower: GAGGAGAAGCCCGGT-TTATCAAGGACACGAAGT

Amplification protocol performed in the Perkin Elmer 2400 thermocycler:

1 time the following step:

pre PCR 2' at 94°C

30 times the following steps:

denaturation 30" at 94°C

annealing 30" at 56°C

elongation 40" at 72°C

1 time the following step:

20 elongation 10' at 72°C

Expected length of specific PCR product: 630 bp.

To quantitatively determine the amount of DNA obtained, amplification products were analysed by electrophoresis on 1% agarose gel in 1xTAE running buffer following standard procedure, as described by Maniatis et al. The samples were compared to a DNA molecular weight marker (MWXVI, Roche).

#### 1.3 Ek/LIC Cloning

The Novagen Ligation Independent Cloning (LIC) kit was used

25

following suppliers instructions in order to obtain directional cloning of the PCR products without the need of restriction enzyme digestion or ligation reactions. The pET-30 Ek/LIC vector, engineered to express the target protein immediately downstream of an enterokinase cleavage site, was chosen.

#### 1.4 Transformation

For good protein expression we chose BL21-Gold(DE3) cells (Stratagene), a derivative of E.coli B, an improved strain of BL21 competent cells. The genotype of the strain is: E.coli B  $F^-$  ompT hsdS ( $r_B^ m_B^-$ ) dcm $^+$  Tet $^-$  gal  $\lambda$ (DE3) endA Hte.

This strain lacks both the lon and the ompT proteases, which can degrade proteins during purification. He phenotype increases the transformation efficiency (>1x108 cfu/µg of pUC18 DNA) In addition the endA gene, that encodes endonuclease I is inactivated (no degradation of plasmid DNA).

In order to obtain competent cells with a high efficiency of transformation, we followed the standard protocol for preparing and electrotransforming BL21-Gold(DE3) cells described in the E. coli Pulser Transformation apparatus manual (BioRad).

Transformation efficiency was tested by using the pUC18 DNA and the pET DNA vectors and the efficiencies obtained were:

- 1x10<sup>10</sup> cfu/μg of pUC18 DNA
- 1x108 cfu/µg of pET DNA

Using these highly electrocompetent cells we were able to obtain a library of approximately 84,000 colonies expressing the randomly mutated Clytin photoprotein.

## 1.5 Plating, induction and charging

Transformed cells were plated on an LB agar plates and grown overnight at 37°C. After overnight colony growth, induction was obtained by

adding 10 mM IPTG and 5 mM EDTA, and incubating for 4 hours at 37°C. Colonies were charged with 10  $\mu$ M Coelenterazine solution and incubated overnight at 4°C in the absence of light.

#### 1.6 CCD camera measurement

5

10

15

20

25

Bioluminescence is assayed by the detection of the signal over a fixed time period of 30" at time 0, after 3 and after 5 minutes from the first measurement.

#### 1.7 Colony picking and re-testing

The best colonies were picked and grown in Iml of LB liquid medium and re-tested using the same experimental conditions described before.

#### 2. In vitro transcription and translation

Translation of the photoproteins was carried out in the rabbit reticulocyte cell-free system (TNT Quick coupled kit, Promega), following the general instructions of the supplier. 500 ng of DNA was used for each *in vitro* transcription/translation reaction mix. The reaction volume (10 µl) included 8 µl of TNT T7 Quick Master Mix, 1.6 µl of DNA, 0.2 µl of the Methionine buffer and 0.2 µl of Coelenterazine (0.5 mM). To this end, 5 µl from each sample of the translation mixture were tested for light emission by the injection of calcium solution and by measuring at the Ascent Luminoskan (Labsystems).

#### 3. Recombinant protein production

For the production of the recombinant protein we followed a small scale protein purification protocol under native conditions. Due to the presence of an N-terminal His tag in our constructs, we decided to purify the expressed proteins on Ni-NTA Spin Columns (Qiagen) following the suppliers' protocol.

#### 4. Calcium concentration curve

In order to evaluate the response of the photoprotein to different calcium concentrations, 0,05 ng/well (96 MTP) of the recombinant protein were charged with Coelenterazine 10 µM overnight at 4°C.

After incubation, different concentrations of CaCl<sub>2</sub> were injected and the light released was measured at the Ascent Luminoskan with an integration time of 20 ms for a total time of 10 seconds.

## 5. Expression of the mutant photoprotein in mammalian cells

#### 5 Reagents

10

25

Restriction enzymes were purchased from New England Biolabs (Beverly, MA) and used according to supplier's instructions. Rapid DNA ligation kit and Fugene transfection reagent were purchased from Roche (Basel, CH). Coelenterazine was from Pharma Tech. International Inc. (Fairfield, NJ). All other chemicals were from standard sources and were of reagent grade or better.

#### Cloning procedure

The most promising mutant photoprotein clones were subcloned for testing their expression in mammalian cells.

2 μl of plasmid were used as template in PCR analysis. In addition a negative control was performed with no template.

Standard PCR procedure were as indicated by Perkin Elmer. PCR protocol was as follows:

#### Primers:

20 Upper primer: TCGTTGGGATCCGCCACCAŢGGCCGACACCGCC

Lower primer: GGGCCCTCTAGATTATCAAGGCACGAA

#### PCR reaction mix:

2 µl	template
5 µl	10 x Pfx Buffer (GIBCO-LifeTechnologies)
1.5 µl	10 mM dNTPs
lμl	50 mM MgSO <sub>4</sub> (GIBCO-LifeTechnologies)
2,5 µl	upper primer (10 μM)
2.5 µl	lower primer (10 μM)

2.5 U Platinum Pfx (GIBCO-LifeTechnologies)

35 μl H<sub>2</sub>O

Amplification protocol performed in Perkin Elmer 2400 thermocycler:

1 time the following step:

5 pre PCR 2' at 94°C

25 times the following steps:

denaturation 15" at 94°C

annealing 30" at 56°C

elongation 40" at 68°C

10 1 time the following step:

15

elongation 10' at 68°C

Expected length of specific PCR product: 630 bp.

Amplification products were analysed by electrophoresis on 1% agarose gel in 1xTAE running buffer following the standard procedure, as described by Maniatis et al.

The PCR product was gel purified using Qiagen columns and digested with BamHI and XbaI restriction enzymes.

pcDNA3neo-/mitoMutated-clytin construction.

An in-house modified pcDNA3 vector (Invitrogen) has been prepared containing the sequence encoding the mitochondrial targeting peptide from subunit VIII of human cytochrome c oxidase (20, 21, 22) so that it could be used in frame at the 5' end of the codon-usage optimized photoprotein gene. The amplification product obtained from the above mentioned PCR has been cloned into this modified pcDNA3 vector lacking the Neomycin resistance gene for expression in mammalian cell lines.

For the mitochondrial targeting the signal sequence is:

 pcDNA3neo-/cytoMutated-clytin construction

The Mutated-clytin gene is obtained from the pcrScript/hMutated-clytin vector by digestion with BamHI and XbaI. The pcDNA3neo- is digested with the BamHI and XbaI restriction enzymes, and purified.

The Mutated-clytin gene is then ligated into the pcDNA3neo-vector to obtain pcDNA3neo-CYTO-hMutated-clytin.

Both constructs obtained were verified by full-length dideoxy sequencing.

#### Cell culture

5

20

25

Culture medium, seeding and incubation: DMEM/F12 with Glutamax (GIBCO cod. 31331-028), 10% FBS, 1% Pen./Strep. (Invitrogen cod. 15140-122), 25 mM Hepes Buffer Solution (GIBCO cod. 15630-056), 1.0 mM Sodium Pyruvate (GIBCO cod. 11360-039), 1.5 g/L Sodium Bicarbonate (GIBCO cod. 25080-060).

Preculture conditions: Cells were seeded for experiments when 70-80% confluent.

Cell culture conditions: Split twice a week: 3.0x10<sup>5</sup> cells/flask T75 (recovery: 8x10<sup>6</sup> cells).

#### Clone Selection Process:

- CHO-K1 were transfected with pcDNA3Neo/MITO-hMutated-clytin.
- 48 h after transfection, the transfected cells were plated in 10x96
   w/p in complete DMEM.
- At confluence the 10x96 w/p were duplicated in 10 white plates. 3 ½ hours before measurements, medium was replaced with 50 μl/well of tyrode (2 mM Ca<sup>2+</sup> and coelenterazine 10 μM).
- Positive clones were selected evaluating:
- First selection was done lysing cells with TRITON X-100. CCD camera conditions: low sensitivity, integration time 1 second reading

for 5 seconds.

5

10

20

25

- Five clones were chosen each diluted in 3x96 w/p.
- At confluence the 15x96 w/p were duplicated in 15 white plates. 3 ½ hours before measurements, medium was replaced with 50 μl/well of tyrode (2 mM Ca<sup>2+</sup> and coelenterazine 10 μM).
- Second selection was done with 10 μM ATP monitoring the kinetics and then cells were lysed with TRITON X-100. CCD camera conditions: low sensitivity, integration time 1 second reading for 30 seconds for the ATP measurement followed by 30 seconds for TRITON X-100.
- Four limiting dilutions of the best clone selected were performed,
   0.3 cells per well in 10x96 w/p.
- The final clone was chosen after the 4th LD selected with ATP
   0.25 μM, coelenterazine 5 μM.
- Complete optimization of the assay was performed on the best responding clone.

#### CCD Camera Measurement parameters:

CHO cells are seeded at different cell densities in 384MTP (500, 750, 1,000, 1,500 cells/well) in growth media supplemented as above and measured with a CCD camera 24 hours and 48 hours after plating. Prior to experiments growth medium is removed and cells are loaded with Tyrode buffer plus coelenterazine at 37°C for 3 hours. Luminescence is finally monitored by CCD camera after the addition of the agonist (30 sec. kinetic).

## Fluorometric Imaging Plate Reader (FLIPR®) Measurements

FLIPR<sup>3</sup> Settings for Assays standard protocol for luminescence detection

- 384 white wall clear bottom plates
- cell plating 24 hrs before the experiment
- · medium removal

- addition of tyrode plus coelenterazine 25 μl/well
- incubation 4 hrs at 37°C
- experiments run at FLIPR<sup>3</sup>: compound (2X) injection in tyrode buffer (25 μl/well).
- All parameter values are the instrumentation default values, except for the following:

#### Pre-assay steps:

1) Camera configuration:

Exposure length = 0.7

10 Camera gain = 200

2) Sequence parameters:

Dispense 384 well head

Height =  $30 \mu l$ 

Speed =  $25 \mu l/sec$ 

15 <u>FLIPR TETRA</u> Settings for Assays standard protocol for luminescence detection

- 384 white wall clear bottom plates
- cell plating 24 hrs before the experiment
- medium removal
- 20 addition of tyrode plus coelenterazine 25 μl/well
  - incubation for 4 hrs at 37°C
  - experiments run at FLIPR<sup>TETRA</sup>: compound (2X) injection in tyrode buffer (25 μl/well).

All parameter values are the instrumentation default values, except for the following:

#### Set up read mode:

1) Camera configuration:

Camera Gain = 200

Exposure time = 0.50

2) Sequence parameters:

Height =  $30 \mu l$ 

Speed =  $25 \mu l/sec$ 

#### 5 Reference Compounds:

15

ATP (Sigma, A-7699) was dissolved in  $H_2O$  at a concentration of 100 mM and stored in aliquots at -20°C. Working solution was freshly prepared in tyrode buffer.

Tyrode Buffer composition: NaCl 130 mM, KCl 5 mM, CaCl<sub>2</sub> 2 mM, 10 MgCl<sub>2</sub> 1 mM, NaHCO<sub>3</sub> 5 mM e HEPES 20 mM, pH 7;4.

### **DESCRIPTION OF THE DRAWINGS**

- Fig. 1: Re-test of mutant colonies with three calcium concentrations.
- Fig. 2: In vitro Transcription & Translation. Measurement of light emission upon 5 pM calcium injection.
- Fig. 3: In vitro Transcription & Translation. Measurement of light emission upon 1 mM calcium injection.
- Fig. 4: Kinetics of the measurement of light emission upon 1 mM calcium injection for the recombinant photoproteins.
- Fig. 5: Calcium dose-response curve of the recombinant photoprotein 20 25N03b.
  - Fig. 6: Peak light intensity from recombinant photoproteins upon 1 mM calcium injection.
  - Fig. 7: Kinetics of 10 mM ATP response at the CCD camera for the CHOK1/mito25N03b cell line.
- Fig. 8: ATP dose-response curve of the CHOK1/mito25N03b cell line.

#### **EXAMPLES**

1. Generation of a Randomly Mutagenized Library and Screening

The Random mutant Library was obtained by using the GeneMorph II

Random Mutagenesis kit (Stratagene). In order to achieve a high mutation rate two different initial amounts of target DNA were used: 0.1 ng and 0.01 ng. A total of 83305 bacterial colonies were tested for luminescence activity. Of these, 1089 were positive and therefore had bioluminescent characteristics. The best colonies were picked, for a total of 289 colonies, and of these, 16 resulted best after a re-test. Finally 8 colonies were chosen and analysed.

Figure 1 shows the results obtained in a three test-point CaCl<sub>2</sub> dose-response curve obtained with the 8 mutants.

#### 2. Photoprotein assay

5

10

15

20

25

5  $\mu$ l of the translation mixture were directly mixed with 95  $\mu$ l of PBS solution in a 96 well plate which was mounted into the Luminometer (Berthold). To trigger photoprotein light emission, a 5 pM CaCl<sub>2</sub> solution was injected into the well and luminescence recorded for 10 seconds.

The results of the in vitro transcription and translation of the DNA of the 8 mutants are shown in Fig. 2.

A new in vitro transcription and translation experiment was carried out with the best responding mutant, 25N03b (sequence ID n°15), and with the Acquorin photoprotein (Fig. 3) in order to have a comparison of the light emission upon the injection of 1mM CaCl<sub>2</sub> solution.

#### 3. Recombinant Photoprotein and Calcium concentration curve

Recombinant photoproteins corresponding to some mutants were produced under native conditions following a small scale purification protocol. Light emission was measured upon 1 mM calcium injection and the corresponding kinetics are shown in Fig. 4.

The recombinant mutant photoprotein, corresponding to clone 25N03b, was better characterized with a complete calcium dose-response curve, which can be seen in Fig. 5.

In another experiment, the recombinant photoprotein 25N03b was

compared to recombinant Aequorin. The light emission recorded upon 1 mM CaCl<sub>2</sub> injection was surprisingly higher in the case of the 25N03b mutant as shown in Fig. 6.

#### 4. Cell-based functional assays

5

10

15

20

25

4.1 The CHOmito25N03b-expressing clone (CHOK1/mito25N03b) has been obtained by transfection of CHO-K1 cells (Materials and Methods). 48 hours after transfection the cells were trypsinized and plated into 10x96 MTP (Microtiter Plates) in complete MEM. At confluence, the 10x96 MTP were duplicated using MATRIX (Hudson, NH, USA) in 10x96 white MTP. 3 hours before measurement the medium was replaced with 50 µl/well of tyrode buffer 2 mM Ca2+ and 10 µM coelenterazine. The clones were selected on the basis of their functional response (luminescent signal) to ATP, which is known to stimulate the CHO endogenous receptor P2Y and to rise the cytoplasmic Ca2+ concentration. At the end of each experiment, cells were lysed by perfusion of a solution containing TritonX-100. The active photoprotein was reconstituted incubating the cells with 10 µM coelenterazine diluted in tyrode buffer containing 2 mM Ca<sup>2+</sup>, in the dark, at 37°C, in a 5% CO<sub>2</sub> atmosphere for 3 hrs. For light emission measurement, cells were lysed in the presence of calcium and the emitted luminescence recorded. The number of photons emitted during the first 30 seconds was integrated by a CCD camera and visualized on the screen. Cells transfected with an empty plasmid or untransfected did not increase photon-emission. To detect changes in calcium concentrations, 10 µM ATP was injected, and the kinetics of the response determined. The curve obtained is shown in Fig. 7.

The final clone was seeded in 384 w/p and tested with increasing ATP concentrations as shown in the dose response curve in Fig. 8.

4.2 The CHOK1/mito25N03b cell line was transfected with a G-protein coupled receptor, the Adenosine A3 receptor, and with a chimeric Gα protein,

10

15

20

25

in order to switch the signal to the PLC/IP pathway. A stable cell line was generated (from now on referred as CHOK1/mito25N03b/A3 clone).

Upon stimulation with its agonist, the A3 receptor induces an increase in intracellular calcium concentration which is measured by mito25N03b luminescence.

4.3 The final clones of the CHO cell lines expressing mito25N03b and the human Adenosine A3 receptor were grown to 80-95% confluence in tissue culture flasks and harvested by trypsinization. Cells were dispensed at different cell densities in 384 w/p in growth medium (DMEM/F12 containing 10% Foetal Bovine Serum) and incubated for 24 h and 48h at 37°C in a humidified incubator at 5% CO<sub>2</sub>. On the day of the experiment, the culture medium was removed and, for luminescence experiments, cells were loaded with 5 μM coelenterazine for 3h, at 37°C, 5% CO<sub>2</sub>. Calcium response was stimulated by addition of different concentrations of ATP to each well. The kinetics of flash luminescence was followed using a CCD camera, which injects reagents and records light emission from each single well.

The agonist and antagonist were diluted in tyrode buffer at different concentrations. Approximately 25  $\mu$ l of these solutions were separately injected into each well and the response measured with the CCD camera instrumentation. The emitted light was immediately recorded at different time intervals.

## 5. Cell-based fluorescence assays at the FLIPR<sup>384</sup>

The CHOK1/mito25N03b/A3 are tested at FLIPR<sup>384</sup> by measuring the fluorescence signal induced by the activation of the transfected receptor.

The cells are incubated with the Calcium assay kit (Molecular Devices Corporation, Sunnyvale, CA, USA).

These fluorescence experiments are carried out to evaluate the luminescence measured with the FLIPR instrumentation.

#### REFERENCES

- 1. Kendall, J.M., and Badminton, M.N. (1998) Aequorea victoria bioluminescence moves into an exciting new era Trends Biotechnology. 16(5):216-24.
- 2. Campbell, A.K., Hallet, R.A., Daw, M.E., Ryall, R.C., Hart and Herring P.J. (1981) Application of the photoprotein obelin to the measurament of free Ca<sup>++</sup> in cells. In Bioluminescence and Chemiluminescence, basic Chemistry and Analytical applications (Edited by M.A. deLuca and W.D. McElroy), pp. 601-607. Academy Press, New York.
- Herring, P.J. (1979) Some features of the bioluminescence of the radiolarian Thalassicola sp. Mar. Biol. 53, 213-216.
  - 4. Shimomura, O., Johnson F.H., and Saiga, Y (1962) Extraction, purification and properties of aequorin, a bioluminescent protein from the luminous hydromedusan, *Aequorea*. J. Cell. Comp. Physiol. 59,223-239.
- 5. Shimomura, O., Johnson F.H., and Saiga, Y (1963) Further data on the bioluminescent protein, aequorin. J. Cell. Comp. Physiol. 62, 1-8.
  - 6. Morin, J.G. and Hastings J.W. (1971) Biochemistry of the bioluminescence of colonial hydroids and other coelenterates. J. Cell. Physiol. 77, 305-311.
- Shimomura, O., Johnson, F.H. and Saiga, Y. (1963) Extraction and properties of halistaurin, a bioluminescent protein from the hydromedusan Halistaura. J. Cell. Physiol. 62, 9-15.
  - 8. Shimomura, O., and Shimomura, A. (1985) Halistaurin, phialidin and modified forms of aequorin as Ca<sup>++</sup> indicator in biological systems. Biochem.
- 25 J. 228,745-749.
  - 9. Levine, L.D., and Ward, W.W. (1982) Isolation and characterization of a photoprotein "phialidin" and a spectrally unique green-fluorescent protein from the bioluminescent jellyfish *Phialidium gregarium*. Comp. Biochem.

Physiol. 72B,77-85.

15

- 10. Morin, J.G. and Hastings (1971) Energy transfer in a bioluminescent system. J. Cell. Physiol. 77, 313-318.
- Campbell, A.K. (1974) Extraction, partial purification and properties of
   obelin the calcium-activated protein from the hydroid *Obelia geniculata*.
   Biochem.J. 143,411-418.
  - 12. Ward, W.W. and Selinger (1974) Extraction and purification of calcium-activated photoprotein from the ctenophores *Mnemiopsis sp.* and *Bern ovata*. Biochemistry 13, 1491-1499.
- 10 13. Ward, W.W. and Seliger H.H. (1974) Properties of mnemiopsin, and berovin, calcium-activated photoproteins from the ctenophores *Mnemiopsis* sp. and *Beroë ovata*. Biochemistry 13, 1500-1510.
  - 14. Johnson, F.H. and Shimomura, O. (1978) Introduction to the bioluminescence of medusae, with special reference to the photoprotein aequorin. Methods Enzymol. 57,271-291.
  - 15. Illarionov B.A., Bondar V.S., Illarionova V.A., Vysotski E.S. Sequence of the cDNA encoding the Ca<sup>++</sup>-activated photoprotein obelin from the hydroid polyp *Obelia longissima*. Gene. 1995 14;153(2):273-4.
  - 16. Blinks, J.R., Weir, W.G., Hess, P. and Prendergast, F.G. (1982).
- 20 Measurement of Ca<sup>++</sup> concentrations in living cells. Prog. Biophys. Mol. Biol. 40,1-114.
  - 17. Markova S.V., Vysotski E.S., Blinks J.R., Burakova L.P., Wang B.C., Lee J., (2002) Obelin from the bioluminescent marine hydroid Obelia geniculata: cloning, expression, and comparison of some properties with those
- 25 of other Ca2+-regulated photoproteins. Biochemistry. 2002 Feb 19;41(7):2227-36.
  - 18. Inouye S., Tsuji F.I. (1993) Cloning and sequence analysis of cDNA for the Ca(2+)-activated photoprotein, clytin. FEBS Lett. Jan 11;315(3):343-6.

- 19. Tsuji F.I., Ohmiya Y., Fagan T.F., Toh H., Inouye S. (1995) Molecular evolution of the Ca(2+)-binding photoproteins of the Hydrozoa. Photochem Photobiol. Oct 62(4):657-61.
- 20. Rizzuto, R., Simpson, A.W.M., Brini, M. and Pozzan, T. (1992) Rapid
   5 changes of mitochondrial Ca2+ revealed by specifically targeted recombinant aequorin. Nature, 358, 325-328.
  - 21. Rizzuto, R., Brini, M., Murgia, M. and Pozzan, T. (1993) Microdomains with high Ca2+ close to IP3-sensitive channels that are sensed by neighbouring mitochondria. Science, 262, 744-747.
- 10 22. Rizzuto, R., Bastianutto, C., Brini, M., Murgia, M. and Pozzan, T. (1994) Mitochondrial Ca2+ homeostasis in intact cells. J. Cell Biol., 126, 1183-1194.

#### **CLAIMS**

5

- 1. Isolated photoprotein containing an amino acid sequence which:
  - a) is able to bind coelenterazine and calcium, producing bioluminescence;
  - b) is identical by at least 90% to SEQ ID NO: 1 (Clytin);
  - c) in sequence alignment with SEQ ID 1 (Clytin), presents one of the following single or multiple substitutions (the residue positions are referred to SEQ ID NO: 1):
- i)  $C_{54} \rightarrow S$ ;
  - ii)  $S_{132} \rightarrow C$ ;
  - iii)  $K_{48}\rightarrow R$ ,  $N_{195}\rightarrow D$ ;
  - iv)  $Q_{68}\rightarrow R$ ,  $A_{90}\rightarrow V$ ,  $T_{184}\rightarrow I$ ;
  - v)  $Y_{82} \rightarrow F$ ,  $K_{110} \rightarrow N$ ,  $F_{125} \rightarrow L$ ,  $S_{149} \rightarrow R$ ;
- 15 vi)  $G_{142} \rightarrow C$ ;
  - vii)  $I_{53} \rightarrow V$ ,  $S_{149} \rightarrow R$ ;
  - viii)  $N_{18} \rightarrow D$ ,  $I_{40} \rightarrow V$ ,  $K_{56} \rightarrow R$ ,

a functional derivative or fragment thereof.

- The photoprotein of claim 1, containing an amino acid sequence
   identical by at least 95% to SEQ ID NO: 1.
  - 3. The photoprotein of claim 2, containing an amino acid sequence identical by at least 98% to SEQ ID NO: 1.
  - 4. The photoprotein of claim 3, containing an amino acid sequence which is selected from the group consisting of SEQ ID NO: 2, 3, 4, 5, 6, 7, 8, 9.
- 25 5. A photoprotein according to claims 1-4, wherein said amino acid sequence is fused to a mitochondrial target sequence.
  - 6. An isolated polynucleotide encoding a photoprotein according to claims 1-5.

- 7. The polynucleotide of claim 6, having the sequence of SEQ ID NO: 10, 11, 12, 13, 14, 15, 16, 17.
- 8. An expression vector containing a polynucleotide according to anyone of claims 6-7.
- 5 9. A prokaryotic or eukaryotic host cell containing the vector of claim 8.
  - 10. A mammalian host cell according to claim 9.
  - 11. A method in vitro for detecting changes in intracellular calcium concentration which comprises:
    - a) providing a cell expressing a photoprotein according to claims 1-5;
- b) contacting the cell with an agent stimulating calcium influx or calcium release from intracellular stores;
  - c) detecting the photoprotein bioluminescence.
  - 12. A method of screening compounds modulating intracellular calcium concentration, which comprises:
    - a) providing a cell expressing a photoprotein of claims 1-5;
      - b) contacting the cell with the candidate compound;
      - c) detecting the bioluminescence of the photoprotein.
  - 13. A method according to claims 12 or 13, which is carried out in a high-throughput format.
- 20 14. A method according to claim 13, which is carried out with a high throughput optical screening apparatus suited for multi-sample analysis.
  - 15. The use of a photoprotein according to claims 1-5 as intracellular calcium indicator.
- 16. The use of a photoprotein according to claim 15 in a cell-based high throughput assay.
  - 17. The use of a photoprotein according to claims 1-5 for the preparation of a diagnostic composition.

#### **ABSTRACT**

## PHOTOPROTEINS WITH ENHANCED BIOLUMINESCENCE AND ASSAYS USING THE SAME

5

The present invention relates to photoproteins with enhanced bioluminescence obtained by mutagenesis of clytin, to their use as intracellular calcium indicators and in cell-based assays.

1 SEQUENCE LISTING

. <110> AXXAM SRL 5 <120> photoproteins with enhanced bioluminescence and assays using the same 10 <130> 1489EUR 15 <160> 22 20 <170> PatentIn version 3.1 25 <210> 1 <211> 196 <212> PRT 30 <213> Clytia gregaria 35 <400> 1 Met Ala Asp Thr Ala Ser Lys Tyr Ala Val Lys Leu Arg Pro Asn Phe 40 Asp Asn Pro Lys Trp Val Asn Arg His Lys Phe Met Phe Asn Phe Leu 25 45 Asp Ile Asn Gly Asp Gly Lys Ile Thr Leu Asp Glu Île Val Ser Lys 40 Ala Ser Asp Asp Ile Cys Ala Lys Leu Gly Ala Thr Pro Glu Gln Thr 50 Lys Arg His Gln Asp Ala Val Glu Ala Phe Phe Lys Lys Ile Gly Met 70. 55 Asp Tyr Gly Lys Glu Val Glu Phe Pro Ala Phe Val Asp Gly Trp Lys 60 Glu Leu Ala Asn Tyr Asp Leu Lys Leu Trp Ser Gln Asn Lys Lys Ser

105

Leu Ile Arg Asp Trp Gly Glu Ala Val Phe Asp Ile Phe Asp Lys Asp 115 120 125

- 5 Gly Ser Gly Ser Ile Ser Leu Asp Glu Trp Lys Ala Tyr Gly Arg Ile 130 135
- Ser Gly Ile Cys Ser Ser Asp Glu Asp Ala Glu Lys Thr Phe Lys His 10 145 150 155 160

Cys Asp Leu Asp Asn Ser Gly Lys Leu Asp Val Asp Glu Met Thr Arg 165 170 175

Gln His Leu Gly Phe Trp Tyr Thr Leu Asp Pro Asn Ala Asp Gly Leu

20
Tyr Gly Asn Phe Val Pro

195

25 <210> 2

<211> 198

<212> PRT

<213> Unknown

35 <220>

<223> Clytin mutant: mutClyKl

<400> 2

- 40
  Met Ala Asp Thr Ala Ser Lys Tyr Ala Val Lys Leu Arg Pro Asn Phe
  1 1 5 10 15
- 45 Asp Asn Pro Lys Trp Val Asn Arg His Lys Phe Met Phe Asn Phe Leu
  20 25 30
- Asp Ile Asn Gly Asp Gly Lys Ile Thr Leu Asp Glu Ile Val Ser Lys 35 40 45
- Ala Ser Asp Asp Ile Ser Ala Lys Leu Gly Ala Thr Pro Glu Gln Thr 50 55 60

Lys Arg His Gln Asp Ala Val Glu Ala Phe Phe Lys Lys Ile Gly Met 65 70 75 80

Asp Tyr Gly Lys Glu Val Glu Phe Pro Ala Phe Val Asp Gly Trp Lys

Glu Leu Ala Asn Tyr Asp Leu Lys Leu Trp Ser Gln Asn Lys Lys Ser 100 105 110

- 5 Leu Ile Arg Asp Trp Gly Glu Ala Val Phe Asp Ile Phe Asp Lys Asp
- Gly Ser Gly Ser Ile Ser Leu Asp Glu Trp Lys Ala Tyr Gly Arg Ile
  10 130 135 140
- Ser Gly Ile Cys Ser Ser Asp Glu Asp Ala Glu Lys Thr Phe Lys His 145 150 155 160

Cys Asp Leu Asp Asn Ser Gly Lys Leu Asp Val Asp Glu Met Thr Arg 165 170 175

- 20
  Gln His Leu Gly Phe Trp Tyr Thr Leu Asp Pro Asn Ala Asp Gly Leu
  180
  185
  190
- 25 Tyr Gly Asn Phe Val Pro 195
- <210> 3

<211> 198

<212> PRT

- 35 <213> Unknown
- <220>

<223> Clytin mutant: mutClyK4

<400> 3

- 45 Met Ala Asp Thr Ala Ser Lys Tyr Ala Val Lys Leu Arg Pro Asn Phe 1 5 10
- Asp Asn Pro Lys Trp Val Asn Arg His Lys Phe Met Phe Asn Phe Leu 20 25 30
- Asp Ile Asn Gly Asp Gly Lys Ile Thr Leu Asp Glu Ile Val Ser Lys 35 40 45
  - Ala Ser Asp Asp Ile Cys Ala Lys Leu Gly Ala Thr Pro Glu Gln Thr
    50 55 60
- 60
  Lys Arg His Gln Asp Ala Val Glu Ala Phe Phe Lys Lys Ile Gly Met
  65
  70
  75
  90

Asp Tyr Gly Lys Glu Val Glu Phe Pro Ala Phe Val Asp Gly Trp Lys 5 Glu Leu Ala Asn Tyr Asp Leu Lys Leu Trp Ser Gln Ash Lys Lys Ser 105 110 Leu Ile Arg Asp Trp Gly Glu Ala Val Phe Asp Ile Phe Asp Lys Asp 10 120 Gly Ser Gly Cys Ile Ser Leu Asp Glu Trp Lys Ala Typ Gly Arg Ile 135 15 Ser Gly Ile Cys Ser Ser Asp Glu Asp Ala Glu Lys Thr Phe Lys His 20 Cys Asp Leu Asp Asn Ser Gly Lys Leu Asp Val Asp Glu Met Thr Arg 25 Gln His Leu Gly Phe Trp Tyr Thr Leu Asp Pro Asn Ala Asp Gly Leu Tyr Gly Asn Phe Val Pro 30 195 <210> 4 35 <211> 198 <212> PRT <213> Unknown 40 <220> 45 <223> Clytin mutant: 1F10 mutant <400> 4 Met Ala Asp Thr Ala Ser Lys Tyr Ala Val Lys Leu Arg Pro Asn Phe 50 Asp Asn Fro Lys Trp Val Asn Arg His Lys Phe Met Phe Asn Phe Leu 55 Asp Ile Asn Gly Asp Gly Lys Ile Thr Leu Asp Glu Ile Val Ser Arg 60 Ala Ser Asp Asp Ile Cys Ala Lys Leu Gly Ala Thr Pro Glu Gln Thr

		<u>, 5</u>													
	Lys Arg 65	His	Gln	Asp	Ala 70	Val	Glu	Ala	Phe	Phe 75	Lys	ГÀз	Ile	Gly	80 80
5	Asp Tyr	Gly	Lys	Glu 85	Val	Glu	Phe	Pro	Ala 90	Phe	Val	Asp	Gly	Trp	Ly
10	Glu Leu	Ala	Asn 100	Tyr	Asp	Leu	Lys	Leu 105	Trp	Ser	Gln	Asn	Lys 110	Гўз	Se:
15	Leu lle	Arg 115	Asp	Trp	Gly	Glu	Ala 120	Val	Phe	Asp	Ile	Phe 125		ГÀг	Asj
	Gly Ser 130	Gly	Ser	Ile	Ser	Leu 135	Asp	Glu	Trṗ	Lys	Ala 140	Туг	Gly	Arg	Il
20	Ser Gly 145	Ile	Суѕ	Ser	Ser 150	Asp	Glu	Asp	Ala	Glu 155	Lys	Thr	Phe	Lys	Ні: 160
25	Cys Asp	Fen	Asp	Asn 165	Ser	Gly	Lys	Leu	Asp 170	Val	Asp	Glu	Met	Thr 175	Arg
30	Gln His	Leu	180 Gly	Phe	Trp	Tyr	Thr	Leu 185	Asp	Pro	Asn	Ala	Asp 190	Gly	Let
35	Tyr Gly	Asp 195	Phe	Val	Pro						· · · · · · · · · · · · · · · · · · ·				
	<210>	5										H			
40	_	19 <b>5</b> PRT													
45	<213> (	Unkno	own												
	<220>														
50	<223>	Clyti	in mu	ıtant	: 1E	17 mu	ıtant	:							
	<400> 5														
55	Met Ala 1	Asp	Thr	Ala 5	Ser	гåг	Туг	Ala	Val 10	Lys	Leu	Arg	Pro	Asn 15	Phe
co.	Asp Asn	Pro	Lys 20	Trp	Val	Asn	Arg	His 25	Lys	Phe	Met	Phe	Asn 30	Phe	Leu
60	Asp Ile	Asn 35	Gly	Asp	Gly	Lys	Ile 40	Thr	Leu	Asp	G <b>l</b> u	Ile 45	Val	Ser	Гуз

Ala Ser Asp Asp Ile Cys Ala Lys Leu Gly Ala Thr Pro Glu Gln Thr 50 55 60 4

- 5 Lys Arg His Arg Asp Ala Val Glu Ala Phe Phe Lys Lys Ile Gly Met 65 70 75 80
- Asp Tyr Gly Lys Glu Val Glu Phe Pro Val Phe Val Asp Gly Trp Lys 90 95

Glu Leu Ala Asn Tyr Asp Leu Lys Leu Trp Ser Gln Ash Lys Lys Ser 100 105 110

Leu Ile Arg Asp Trp Gly Glu Ala Val Phe Asp Ile Phe Asp Lys Asp 115

20
Gly Ser Gly Ser Ile Ser Leu Asp Glu Trp Lys Ala Tyr Gly Arg Ile
130
135
140

- Ser Gly Ile Cys Ser Ser Asp Glu Asp Ala Glu Lys Thr Phe Lys His 145 150 155 160
- Cys Asp Leu Asp Asn Ser Gly Lys Leu Asp Val Asp Glu Met Thr Arg 165 170 175

Gln His Leu Cly Phe Trp Tyr Ile Leu Asp Pro Asn Ala Asp Gly Leu 180 185 190

Tyr Gly Asn Phe Val Pro 195

40 <210> 6

<211> 198

45 <212> PRT

<213> Unknown

**50** <220>

<223> Clytin mutant: 1C12 mutant

55 <400> 6

Met Ala Asp Thr Ala Ser Lys Tyr Ala Val Lys Leu Arg Pro Asn Phe 1 5 10 15

Asp Asn Pro Lys Trp Val Asn Arg His Lys Phe Met Phe Asn Phe Leu 20 25 30

Asp Ile Asn Gly Asp Gly Lys Ile Thr Leu Asp Glu Ile Val Ser Lys 35 40 ; 45

- 5 Ala Ser Asp Asp Ile Cys Ala Lys Leu Gly Ala Thr Pro Glu Gln Thr 50 55 60
- Lys Arg His Gln Asp Ala Val Glu Ala Phe Phe Lys Lys Ile Gly Met 65 70 75  $_{\ast}$  80
- Asp Fhe Gly Lys Glu Val Glu Phe Pro Ala Phe Val Asp Gly Trp Lys 85 90 95
  - Glu Leu Ala Asn Tyr Asp Leu Lys Leu Trp Ser Gln Asn Asn Lys Ser 100 105 110
- Leu Ile Arg Asp Trp Gly Glu Ala Val Phe Asp Ile Leu Asp Lys Asp 115 120 125
- 25 Gly Ser Gly Ser Ile Ser Leu Asp Glu Trp Lys Ala Tyr Gly Arg Ile 130 135 140
- Ser Gly Ile Cys Arg Ser Asp Glu Asp Ala Glu Lys Thr Phe Lys His 150 155 160
- Cys Asp Leu Asp Asn Ser Gly Lys Leu Asp Val Asp Glu Met Thr Arg 165 170 175
  - Gln His Leu Gly Phe Trp Tyr Thr Leu Asp Pro Asn Ala Asp Gly Leu 180 185 4 190
- 40
  Tyr Gly Asn Phe Val Pro
  195
- 45 <210> 7
  - <211> 198
- <212> PRT
  - <213> Unknown
- 55 <220>
  - <223> Clytin mutant: 25N03b mutant
- <400> 7
  - Met Ala Asp Thr Ala Ser Lys Tyr Ala Val Lys Leu Arg Pro Asn Phe 1 5 10 15

ጵ

Asp Asn Pro Lys Trp Val Asn Arg His Lys Phe Met Phe Asn Phe Leu 20 25 30

- 5 Asp Ile Asn Gly Asp Gly Lys Ile Thr Leu Asp Glu Ile Val Ser Lys 35 40 45
- Ala Ser Asp Asp Ile Cys Ala Lys Leu Gly Ala Thr Pro Glu Gln Thr  $10 \ 50 \ 55 \ 60$
- Lys Arg His Gln Asp Ala Val Glu Ala Phe Phe Lys Lys Ile Gly Met 65 70 75 80

Asp Tyr Gly Lys Glu Val Glu Phe Pro Ala Phe Val Asp Gly Trp Lys

- Glu Leu Ala Asn Tyr Asp Leu Lys Leu Trp Ser Gln Asn Lys Lys Ser 100 105 110
- 25 Leu Ile Arg Asp Trp Gly Glu Ala Val Phe Asp Ile Phe Asp Lys Asp 115 120 125
- Gly Ser Gly Ser Ile Ser Leu Asp Glu Trp Lys Ala Tyr Cys Arg Ile 30 130 135 140
- Ser Gly Ile Cys Ser Ser Asp Glu Asp Ala Glu Lys Thr Phe Lys His 145 150 155 160
  - Cys Asp Leu Asp Asn Ser Gly Lys Leu Asp Val Asp Glu Met Thr Arg 165 170 175
- Gln His Leu Gly Phe Trp Tyr Thr Leu Asp Pro Asn Ala Asp Gly Leu
  180 185 190
- 45 Tyr Gly Asn Phe Val Pro 195
- <210> B

<211> 198

<212> PRT

- 55 <213> Unknown
- <220>
  60
  <223> 'Clytin mutant: 3Cl2 mutant

g

<400> 8

5

Met Ala Asp Thr Ala Ser Lys Tyr Ala Val Lys Leu Arg Pro Asn Phe 1 5 10 15

Asp Asn Pro Lys Trp Val Asn Arg His Lys Phe Met Phe Asn Phe Leu 20 25 30

Asp fle Asn Gly Asp Gly Lys Ile Thr Leu Asp Glu Ile Val Ser Lys
35 40 45

15 Ala Ser Asp Asp Val Cys Ala Lys Leu Gly Ala Thr Pro Glu Gln Thr 50 55 60

Lys Arg His Gln Asp Ala Val Glu Ala Phe Phe Lys Lys Ile Gly Met 20 65 70 75 80

Asp Tyr Gly Lys Glu Val Glu Phe Pro Ala Phe Val Asp Gly Trp Lys
85 90 95

Glu Leu Ala Asn Tyr Asp Leu Lys Leu Trp Ser Gln Asn Lys Lys Ser 100 105 110

Leu Ile Arg Asp Trp Gly Glu Ala Val Phe Asp Ile Phe Asp Lys Asp 115 120 125

35 Gly Ser Gly Ser Ile Ser Leu Asp Glu Trp Lys Ala Tyr Gly Arg Ile 130 135 140

Ser Gly Ile Cys Arg Ser Asp Glu Asp Ala Glu Lys Thr Phe Lys His 40 145 150 150 155

Cys Asp Leu Asp Asn Ser Gly Lys Leu Asp Val Asp Glu Met Thr Arg 165 170 175

Gln His Leu Gly Phe Trp Tyr Thr Leu Asp Pro Asn Ala Asp Gly Leu 180 185 190

50 Tyr Gly Asn Phe Val Pro 195

55 <210> 9

45

<211> 198

<212> PRT

<213> Unknown

<220>

<223> Clytin mutant: 6H22 mutant

- 5 <400> 9
  - Met Ala Asp Thr Ala Ser Lys Tyr Ala Val Lys Leu Arg Pro Asn Phe 1 5 10 15
- Asp Asp Pro Lys Trp Val Ash Arg His Lys Phe Met Phe Ash Phe Leu 20 25 30
- 15 Asp Ile Asn Gly Asp Gly Lys Val Thr Leu Asp Glu Ile Val Ser Lys 35 40
- Ala Ser Asp Asp Ile Cys Ala Arg Leu Gly Ala Thr Pro Glu Gln Thr 20 50 55 60
- Lys Arg His Gln Asp Ala Val Glu Ala Phe Phe Lys Lys Ile Gly Met 65 70 75 80
  - Asp Tyr Gly Lys Glu Val Glu Phe Pro Ala Phe Val Asp Gly Trp Lys 85 90 95
- Glu Leu Ala Asn Tyr Asp Leu Lys Leu Trp Ser Gln Asn Lys Lys Ser 100 105 110
- 35 Leu Ile Arg Asp Trp Gly Glu Ala Val Phe Asp Ile Phe Asp Lys Asp 115 . 120 . 125
- Gly Ser Gly Ser Ile Ser Leu Asp Glu Trp Lys Ala Tyr Gly Arg Ile
  40 130 135 140
- Ser Gly Ile Cys Ser Ser Asp Glu Asp Ala Glu Lys Thr Phe Lys His 145 150 150 155
  - Cys Asp Leu Asp Asn Ser Gly Lys Leu Asp Val Asp Glu Met Thr Arg
- Gln His Leu Gly Phe Trp Tyr Thr Leu Asp Pro Asn Ala Asp Gly Leu 180 185 190
- 55 Tyr Gly Asn Phe Val Pro 195
- 60 <210> 10 <211> 600

11

<212> DNA <213> Unknown 5 <220> <223> Clytin mutant: mutClyK1\_dna 10 <400> 10 atggccgaca ccgccagcaa gtacgccgtg aagctgaggc ccaacttcga caaccccaag 120 tgggtgaacc ggcacaagtt catgttcaac ttcctggaca tcaacggcga cggcaagatc 15 accotggacg agategtgag caaggecage gacgacatet gegecaaget gggegecaec 180 cccgagcaga ccaagagaca ccaggacgcc gtggaggcct tettcaagaa gatcggcatg 240 20 300 gactacggca aggaggtgga gttccccgcc ttcgtggacg gctggaagga gctggccaac taccacctga agctgtggag ccagaacaag aagagcctca tcagggactg gggcgaggcc 360 gtgttcgaca tettegacaa ggacggeage ggetgeatca geetggatga gtggaaggee 420 25 tacggcagaa tcagcggcat ctgcagcagc gacgaggacg ccgaaaaagac cttcaagcac 480 tgcgacctgg acaacagcgg caagetggac gtggacgaga tgaccagaca gcacctggac 540 30 ttetggtaca ccetggacce caatgeegae ggcctgtacg geaacttegt geettgataa 600 <210> 11 35 <211> 600 <212> DNA <213> Unknown 40 <220> 45 <223> Clytin mutant: mutClyK4\_dna atggeogaca cogecagoaa gtacgoogtg aagotgaggo coaacttoga caaccooaag 60 50 tgggtgaacc ggcacaagtt catgttcaac ttcctggaca tcaacggcga cggcaagatc 1.20 accetggacg agategtgag caaggecage gacgacatet gegetaaget gggegecace 160 cccgagcaga ccaagagaca ccaggacgcc gtggaggcct tcttbaagaa gatcggcatg 240 55 gactacggca aggaggtgga gttccccgcc ttcgtggacg gctgggagga gctggccaac 300 tacgacctga agctgtggag ccagaacaag aagagcctca tcagggactg gggcgaggcc 360 60 gtgttcgaca tcttcgacaa ggacggcago ggctgcatca gcctggatga gtggaaggcc 120 tacggcagaa tcagcggcat ctgcagcagc gacgaggacg ccgaaaaaac cttcaagcac 480

	tgcgacc	tgg	acaacagegg	caagctggac	gtggacgaga	tgaccagaca	geacetggge	540
	ttctggt	aca	ccctggaccc	caatgeegae	ggeetgtacg	gcaacttcgt	gccttgataa	600
5	<210> <211>	12 500						
10								•
10	<212>	DNA		• •	•			÷
	<213>	Unkı	JOWN	•				
15	<220>							
	<223>	Cly	tin mutant:	1F10 mutan	t_dna			
20	<400> atggcc	12 gaca	ccgccagcaa	gtacgccgtg	aagetgagge	ccaactccga	caaccccaag	60
	tgggtg	aacc	ggcacaagtt	catgttcaac	ttectggaca	tcaacggcga	eggeaagate	120
25	accetg	gacg	agatcgtgag	cagggccagc	gacgacatet	gcgccaagct	gggcgcca <b>c</b> c	180
	cccgag	caga	ccaagagaca	ccaggacgcc	gtggaggcct	tottoaagaa	gatcggcatg	240
30	gactac	ggca	aggaggtgga	gttccccgcc	ttcgtggacg	gctggaagga	gctggccaac	300
	tacgac	ctga	agctgtggag	ccagaacaag	aagageetca	tcagggactg	åddedaddee	360
	gtgttc	gaca	tcttcgacaa	ggacggcagc	ggcagcatca	gcctggatga	gtggaaggcc	420
35	tacggc	agaa	tcagcggcat	ctgcagcago	gacgaggacg	ccgaaaagac	cttcaagcac	480
	tgcgac	ctgg	acaacagcgg	caagctggac	gtggacgaga	. tgacçagaca	geacctgggc	540
40	ttctgg	taca	ccctggaccc	caacgecgac	: ggcctgtacg	gegaettegt	gccttgataa	600
	<210>	13		•				
<b>.</b>	<211>	600						
45	<212>	DNA						
	<213>	Unk	nown					
50								
	<2,20>							
55	<223>	Cly	tin mutant:	1H7 mutant	_dna			
	<400> atggcc	13 gaca	ccgccagcaa	gtacgccgtç	, aagctgaggo	ccaacttcga	caaccccaag	60
60	tgggtg	aacc	ggcacaagtt	catgttcaat	: tteatggaca	. tcaacggcga	eggcaagate	120
	accctg	gacg	agatcgtgag	caaggccago	gacgacatct	gegeeaaget ,	gggcgccacc	180

1	3
~	~

				10			
	cccgagcaga	ccaagagaca	ccdddacdcc	gtggaggcct	tcttcaagaa	gatcggcatg	240
	gactacggca	aggaggtgga	gttccccgtc	ttcgtggacg	gctggaagga	gctggccaac	300
5	tacgacctga	agctgtggag	ccagaacaag	aagagcctca	tcagggactg	gggcgaggcc	360
	gtgtttgaca	tettegacaa	ggacggcagc	ggcagcatta	gcctggatga	gtggaaggcc	420
10	tacggtagaa	tcagcggcat	ctgcagcagc	gacgaggacg	ccgaaaagac	cttcaagcac	480
	tgcgacctgg	acaacagcgg	caagctgġac	gtggacgaga	tgaccagaca	gcacctggge	540
	ttctggtaca	teetggaeee	caacgccgac	ggcctgtacg	gcaacttcgt	gccttgataa	600
15	<210> 14				•		
	<211> 600						
20	<211> 000 <212> DNA						
20	<213> Unkn	own					
	(ZIJ) OIIMI						
25	<220>						•
		in mutant:	1C12 mutant	t dina			
30	<400> 14						
	atggccgaca	ccgccagcaa	gtacgccgtg	aagctgaggc	ccaacttcga	caaccccaag	60
	tgggtgaacc	ggcacaagtt	catgttcaac	tteetggaca	tcaacggcga	cggcaagatc	120
35	accctggacg	agatcgtgag	caaggccagc	gacgacatct	gcgccaagct	gggcgccacc	180
	cccgagcaga	ccaagagaca	ccaggacgee	gtggaggcct	tcttcaagee	gateggeatg	240
40	gacttcggca	aggaggtgga	gttccccgcc	ttcgtggacg	gctggaagga	gctggccaac	300
	tacgacctga	agctgtggag	ccagaacaat	aagagcetea	tcagggactg	gggcgaggcc	360
	gtgttcgaca	teetegacaa	ggacggcagc	ggcagcatca	gcctggatga	gtggaaggcc	420
45	tacggcagaa	tcagcggcat	ctgcagaagc	gacgaggacg	ccgaaaagac	cttcaagcac	480
	tgcgacctgg	acaacagcgg	caagctggac	gtggacgaga	tgaccagaca	gcacctgggc	540
50	ttctggtaca	ccctggaccc	caacgecgac	ggcctgtacg	gcaacttcgt	gccttgataa	600
	<210> 15						
55	<211> 600				3	,	
	<212> DNA						
	<213> Unkn	own					
60	•						

	, 14					
	<220>					
	<223> Clytin mutant: 25N03b mutant_dna					
. 5	<400> 15 atggccgaca ccgccagcaa gtacgccgtg aagctgaggc ccaacttcga caaccccaag	60				
10.	tgggtgaace ggcacaagtt catgttcaac ttcctggaca tcaacggcga cggcaagatc	120				
	accetggacg agategtgag caaggeeage gaegacatet gegeeaage gggegeeaee	180				
	cccgagcaga ccaagagaca ccaggacgcc gtggaggcct tottcaagaa gatcggcatg	240				
15	gactacggca aggaggtgga gttccccqcc ttcgtggacg gctggaagga gctggccaac	300				
	tacgaeetga agetgtggag ccagaacdag aagageetea tcagggaetg gggegaggee	360				
	gtgttcgaca tettcgacaa ggacggcagc ggcagcatca gcctggatga gtggaaggcc	420				
20	tactgcagaa tcagcggcat ctgcagcagc gacgaggacg ccgaaaagac cttcaagcac	480				
	tgcgacctgg acaacagcgg caagctggac gtggacgaga tgaccagaca gcacctgggc	540				
25	ttetggtaca ecetggacce caacgeegae ggeetgtaeg geaacttegt geettgataa	600				
	<210> 16					
30	<211> 600					
	<212> DNA					
	<z13> Unknown</z13>					
35						
	<220>					
40	<223> Clytin mutant: 3C12 mutant_dna					
	<400> 16 atggccgaca ccgccagcaa gtacgccgtg aagctgaggc ccaacttcga caaccccaag	60				
45	tgggtgaacc ggcacaagtt catgttcaac ttcctggaca tcaacggcga cggcaagatc	120				
	accotggacg agatogtgag caaggocage gacgacgtot gegecaaget gggegecaec	180				
50	cccgagcaga ccaagagaca ccaggacgcc gtggaggeet tettcaagaa gatcggcatg	240				
	gactacggca aggaggtgga gttccccgcc ttcgtggacg gctggaagga gctggccaac	300				
	tacgacctga agctgtggag ccaaaacaag aagagcctca tcagggactg gggcgaggec	360				
55	gtgttcgaca tcttcgacaa ggacggcagc ggcagcatca gcctggacga gtggaaggcc	420				
	tacggcagaa tcagcggcat ctgcagaagc gacgaggacg ccgaaaagac cttcaagcac	480				

tgcgacctgg acaacagcgg caagctggac gtggacgaga tgaccagaca gcacctgggc

ttctggtaca ccctggaccc caacgccgac ggcctgtacg gcaacttcgt gccttgataa

<211> 5 <212> DNA

<400> 17

<210>

<213> Unknown

600

10 <220>

15

<223> Clytin mutant: 6H22 mutant\_dna

atggccgaca ccgccagcaa gtacgccgtg aagctgaggc ccaacttcga cgaccccaag 60 tgggtgaacc ggcacaagtt catgttcaac ttcctggaca tcaacggcga cggcaaggte 120 20 accetggacg agategtgag caaggecage gacgacatet gegecagget gggegecace 180 ceegageaga ceaagagaea eeaggaegee gtggaggeet tetteaagaa gateggeatg 240 gactacygca aagaggtgga gttccccgcc ttcgtggacg gctggaagga gctggccaac 300 25 tacgacctga agctgtggag ccagaacaag aagagcctca tcagggactg gggcgaggcc 360 gtgttcgaca tettcgacaa ggacggcagc ggcagcatca gcctggatga gtggaaggcc 420 30 tacggcagaa tcagcggcat ctgcagcagc gacgaggacg ccgaaaagac cttcaagcac 480 tgcgacctgg acaacagcgg caagctggac gtggacgaga tgaccagaca gcacctgggc 540

ttctggtaca ccctggaccc caacgccgac ggcctgtacg gcaacttcgt gccttgataa

15

<Z10> 18

<211> 32 40

<212> DNA

<213> Unknown

45

35

<220>

<223> synthetic primer 50 <400> 18

gatgacgacg acaagatggc cgacaccgcc ag

55 <210> 19

<211> 33


<212> DNA

60 <213> Unknown 32

			16	·
	<220>			
	<223>	synthetic primer	e e	
5	<400> gaggag	19 gaage eeggtttate aaggacaega agt	**************************************	. 33
10	<210>	20		
	<211>			
	<212>	DNA		
15	<213>	Unknown		
20	<220>		*	•
	<223>	synthetic primer	-	
25	<400> tegttg	20 ggat cogccaccat ggccgacacc gcc		
	<210>	21		
30	<211>	27		
30	<212>	DNA		
	<213>	Unknown		
35				
	<220>			
40	<223>	synthetic primer		
40	<400> gggccc	21 tota gattatoaag goacgaa		27
45	<210>	22		
	<211>	99		
50	<212>	DNA		
	<213>	Homo sapiens		
55	<400> atgtecq	•	gacag geteggeeeg	geggeteeca 60
	gtgccgd	egog ccaagatoca ttogttggga toogo	cacc	99
60				


Sheet 1/8

FIG. 1

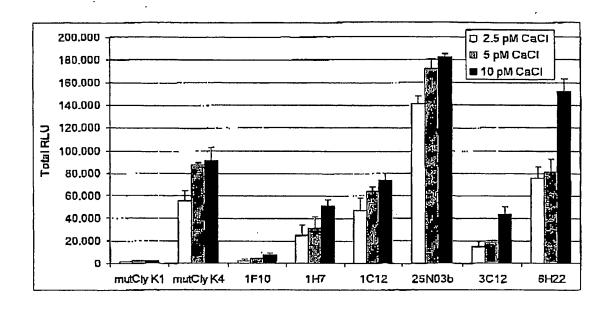


Sheet 2/8

FIG. 2

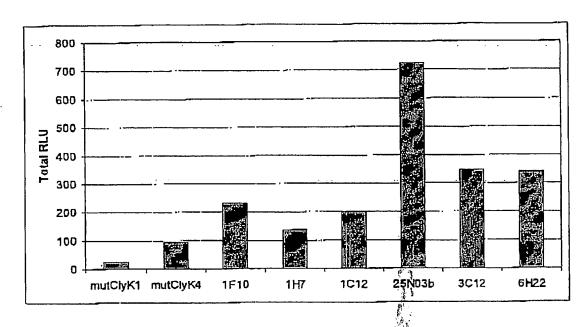


Sheet 3/8



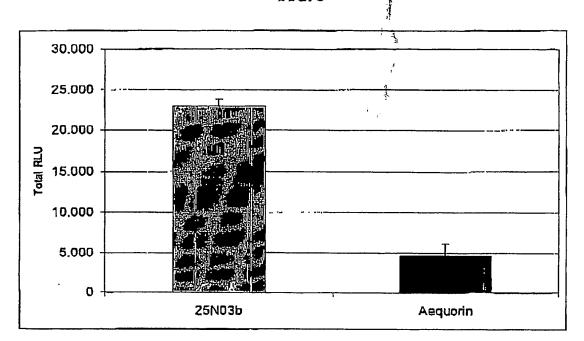


Sheet 4/8

FIG. 4

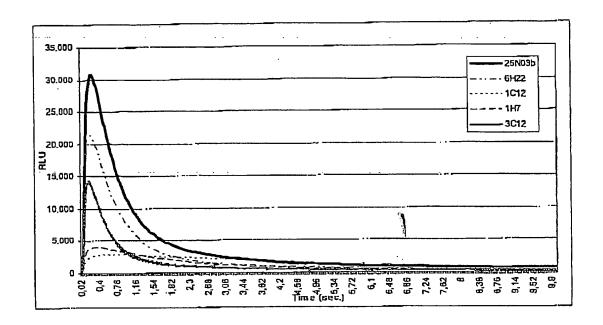


Sheet 5/8

**FIG.** 5

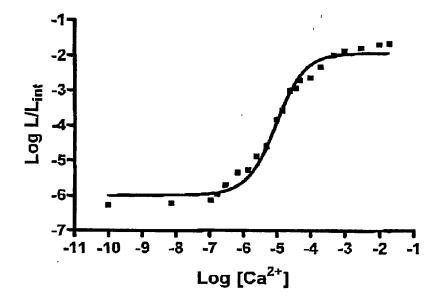


Sheet 6/8

FIG. 6

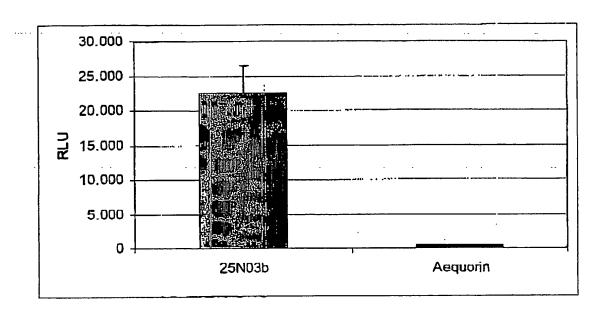


Sheet 7/8

**FIG.** 7

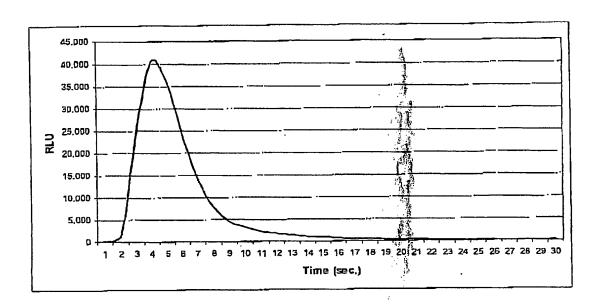

Sheet 8/8

FIG. 8

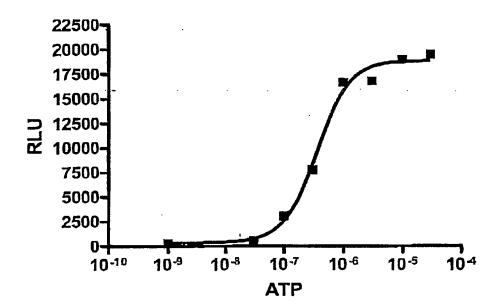

## Document made available under the Patent Cooperation Treaty (PCT)

International application number: PCT/EP2006/002172

International filing date:

09 March 2006 (09.03.2006)

Document type:

Certified copy of priority document

Document details:

Country/Office: EP

Number:

06000171.6

Filing date:

05 January 2006 (05.01.2006)

Date of receipt at the International Bureau:

20 June 2006 (20.06.2006)

Remark:

Priority document submitted or transmitted to the International Bureau in

compliance with Rule 17.1(a) or (b)



## **Best Available Copy**



Europäisches Patentamt GD1

European Patent Office DG1

Office européen des brevets

PCT/EP200 6/002172

09.03.06

Bescheinigung

Certificate

Attestation

Die angehefteten Unterlagen stimmen mit der ursprünglich eingereichten Fassung der auf dem nächsten Blatt bezeichneten europäischen Patentanmeldung überein. The attached documents are exact copies of the European patent application described on the following page, as originally filed.

Les documents fixés à cette attestation sont conformes à la version initialement déposée de la demande de brevet européen spécificée à la page suivante.

Patentanmeldung Nr.

Patent application No.

Demande de brevet n°

06000171.6 / EP06000171

The organization code and number of your priority application, to be used for filing abroad under the Paris Convention, is EP06000171

Der Präsident des Europäischen Patentamts; Im Auftrag

For the President of the European Patent Office Le President de l'Office européen des brevets p.o.

R.C. van Dijk

## **Best Available Copy**



Europäisches Patentamt GD1

**European Patent Office** DG<sub>1</sub>

Office européen des brevets DG<sub>1</sub>

PCT/EP200 6/002172

09.03.06

Anmeldung Nr:

Application no.: Demande no:

06000171.6

Anmeldetag: Date of filing: Date de dépôt:

05.01.06

Anmelder/Applicant(s)/Demandeur(s):

AXXAM S.r.I. Via Massena, 12/7 20145 Milano/IT

Bezeichnung der Erfindung/Title of the invention/Titre de l'invention: (Falls die Bezeichnung der Erfindung nicht angegeben ist, siehe Beschreibung. If no title is shown please refer to the description. Si aucun titre n'est indiqué se referer à la déscription.)

Photoproteins with enhanced bioluminescence and assays using the same

In anspruch genommene Priorität(en) / Priority(ies) claimed / Priorité(s) revendiquée(s) Staat/Tag/Aktenzeichen / State/Date/File no. / Pays/Date/Numéro de dépôt:

Internationale Patentklassifikation / International Patent Classification / Classification internationale de brevets:

C07K14/435

Am Anmeldetag benannte Vertragstaaten / Contracting states designated at date of filing / Etats contractants désignées lors du dépôt:

AT BE BG CH CY CZ DE DK EE ES FI FR GB GR HU IE IS IT LI LT LU LV MC NL PL PT RO SE SI SK TR